

INFUSION SYSTEMS


PRODUCT GUIDE

TIVA Sets & A Sets



TIVA Sets


Total intravenous anaesthesia from pump to patient

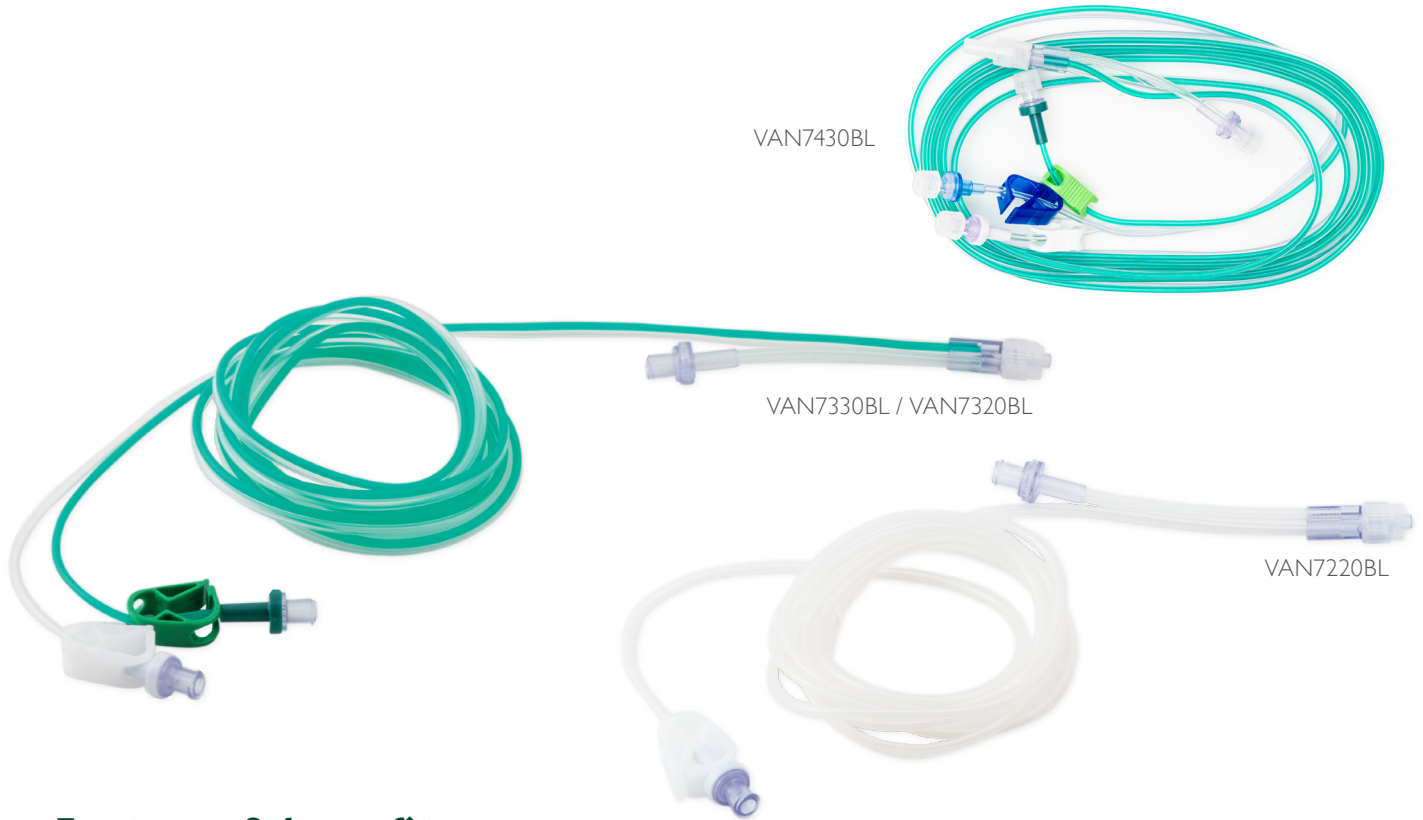

Features & benefits



Colour-coded Siamese tubing
For clearer identification of high and low flow drug administration.

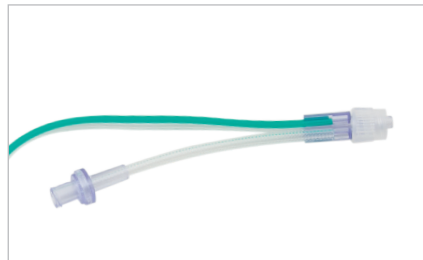

Multi-lumen from hub to tip
Keeps infusions separate, reducing the risk of drugs mixing.




Armoured tubing
Eliminates the risk of kinking, ensuring continued infusions.

Ordering codes

Vygon codes	NHSSC codes	Number of lumen	Gravity line		High flow drug line		Low flow drug line		Low flow drug line		Box qty
			Length (cm)	Priming vol (mL)	Length (cm)	Priming vol (mL)	Length (cm)	Priming vol (mL)	Length (cm)	Priming vol (mL)	
VAN7430BL	-	Quad	8.0	0.5	300	7.0	300	2.5	300	2.5	50
VAN7420BL	-	Quad	8.0	0.5	200	4.8	200	1.8	200	1.8	50
VAN7320BL	FSB1337	Triple	8.0	0.5	200	4.8	200	1.8	-	-	50
VAN7330BL	FSB1338	Triple	8.0	0.5	300	7.0	300	2.5	-	-	50
VAN7230BL	-	Double	8.0	0.5	300	7.0	-	-	-	-	50
VAN7220BL	FSB1336	Double	8.0	0.5	200	4.8	-	-	-	-	50

A Sets

Total protection from pump to patient

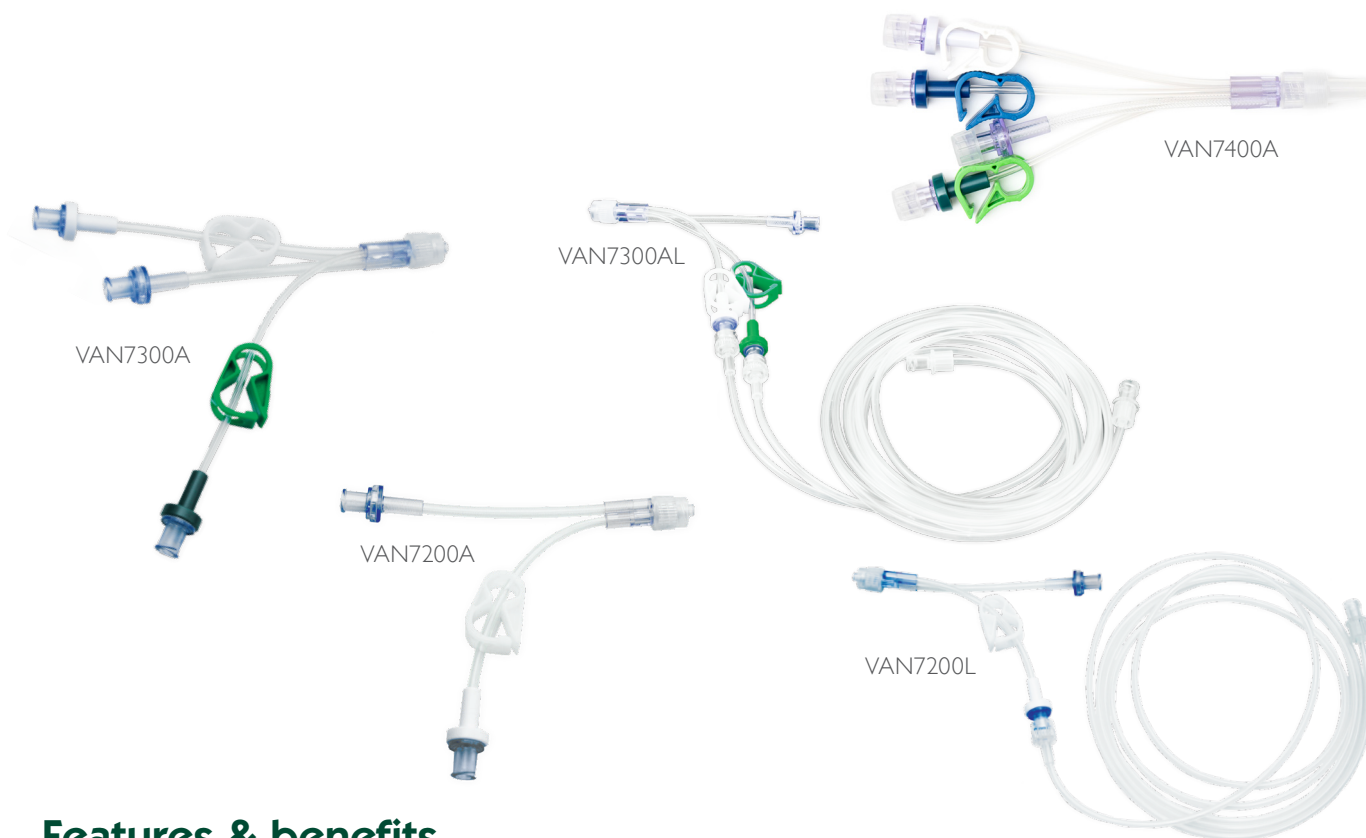

Features & benefits



Anti-syphon valve

One-way valves that require a far greater driving pressure to open, minimising the risk of downloading.



Anti-reflux valve

Minimises the risk of back-tracking and enhances patient safety.



Halkey Roberts clamp

Secure, easy-to-use clamps, which enable one-handed operation during medication changes.

Ordering codes

Vygon codes	NHSSC codes	Number of lumen	Detachable extension line	Gravity line		High flow drug line		Low flow drug line		Low flow drug line		Box qty
				Length (cm)	Priming vol (mL)	Length (cm)	Priming vol (mL)	Length (cm)	Priming vol (mL)	Length (cm)	Priming vol (mL)	
VAN7400A	-	Quad	No	8.0	0.6	10	0.6	10	0.3	10	0.3	50
VAN7300A	FSB1081	Triple	No	8.0	0.5	10	0.6	12	0.3	-	-	50
VAN7300AL	FSB1082	Triple	Yes	8.0	0.5	185	5.0	187	4.7	-	-	50
VAN7200A	FSB1079	Double	No	8.0	0.5	10	0.6	-	-	-	-	50
VAN7200L	FSB1080	Double	Yes	8.0	0.5	210	5.1	-	-	-	-	50



08/2024 - DX|B0005002 V00566 v2

FOR FURTHER INFORMATION, PLEASE CONTACT: info@vygon.co.uk

The specifications shown in this leaflet are for information only and are not, under any circumstances, of a contractual nature. This brochure has been printed on responsibly sourced and sustainable material. To help us reduce our carbon footprint all of our literature is available electronically either from your Product Specialist, on our website or by emailing info@vygon.co.uk

VYGON (UK) LTD, THE PIERRE SIMONET BUILDING, V PARK,
GATEWAY NORTH, LATHAM ROAD, SWINDON, WILTSHIRE SN25 4DL
RECEPTION: +44 (0)1793 748800 WWW.VYGON.CO.UK

 [@vygonuk](https://www.instagram.com/vygonuk)

