

ONCOLOGY

RANGE



The high security connection
for systemic anti-cancer
therapy



Value Life



Stay safe with Qimono

It is widely documented that exposure to Systemic Anti-Cancer Therapy (SACT), can pose a significant health risk to patients during treatment as well as healthcare workers involved in the preparation and administration of the drugs¹.

Many SACT agents are known to be carcinogenic, teratogenic and mutagenic.


Exposure may be through skin contact, skin absorption, inhalation of aerosols and drug particles as well as ingestion and needlestick injuries. This can occur at various points of the treatment process: during the preparation of the medication, whilst the drugs are administered and after therapy².

Adverse health effects reported by chemotherapy nurses include headaches, dizziness, nausea and hair loss³.

Closed systems to minimise the risk

In the COSHH Regulations 2002, cytotoxic drugs are defined as hazardous substances.

The 2015 guidance issued by the Health and Safety Executive (HSE) to minimise the risk to clinicians prioritises the use of totally enclosed systems where reasonably practicable².

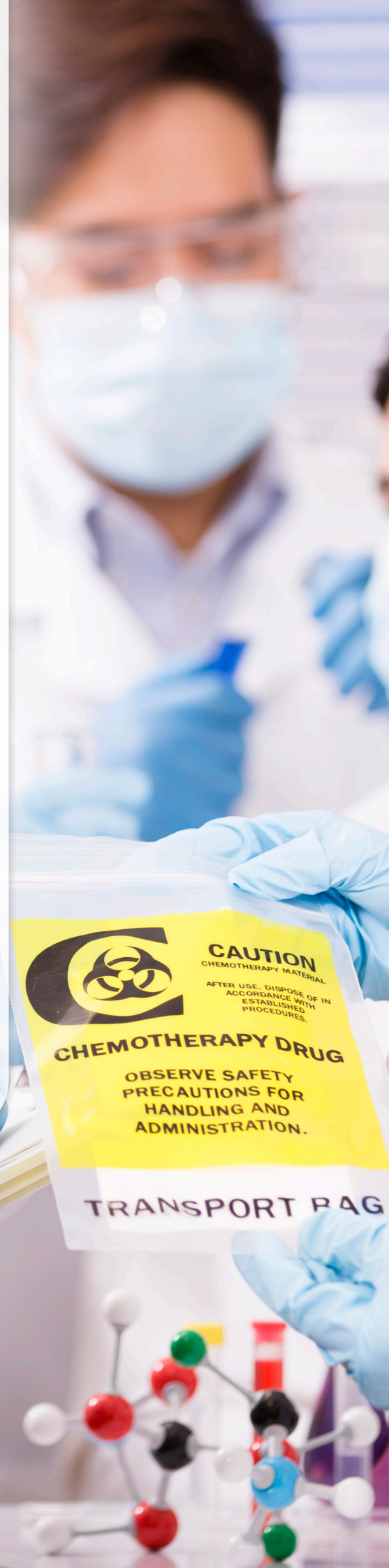

The seven key benefits of Qimono

Qimono is a totally closed system drug transfer device specifically designed to be safe and easy to use by pharmacists and nurses.

The seven key benefits ensure there are secure connections in place for the preparation, transfer and infusion of medication to protect clinicians throughout the delivery of chemotherapy to patients.



1 Luer bayonet provides connections without screwing

- Eliminates the problems linked to the clamping force applied on medical devices
- Reduces the risk of repetitive strain injury for users
- Removes the likelihood of over, or under tightening.

2 Connectable to standard Luer devices

- Connects easily to standard devices such as syringes and extension lines.

3 The two tube design of the internal mechanism provides a direct fluid pathway

- Easy withdrawals and injections
- Allows effective flushing
- Limits ward pump alarms.

4 Chemical resistant materials

- Avoids medical devices cracking.

5 Valves automatically close on disconnection ensuring a barrier against toxic drugs and bacteria

- Protects users from chemical contamination
- Shields drugs from bacterial contamination.

6 Flat membranes mean the easy disinfection of both male ♂ and female ♀ hubs



- Makes it simple to disconnect and connect devices when needed.

7 A complete range of devices for the preparation, transfer and infusion of medication focuses on safety from start to finish




- Makes all handling of cytotoxic drugs safe
- Guarantees that all devices are compatible.

PREPARATION DEVICES FOR MEDICATIONS/SOLUTIONS


Bag spikes

Vygon code	Image	Description	Indication	Quantity box/case
00721801		Qimospike B is an IV fluid bag access spike fitted with Qimo female ♀	Aspiration or addition of medication/fluid from/into an IV fluid bag	50/600
00721811		Qimospike B valve is an IV fluid bag access spike fitted with Qimo female ♀ and non-return valve	Aspiration only from an IV fluid bag, thus avoiding the accidental contamination of the IV solution	25/450

Vial

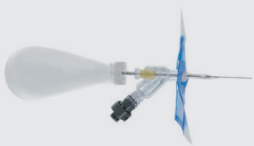

Vygon code	Size	Image	Description	Indication	Quantity box/case
00721901	ø14		Qimocap is a vial access cap fitted with Qimo female ♀. It has been specifically designed to create a permanent connection with the medication vial	Withdrawal and reconstitution of small volumes of medication in completely closed environment	50/900
00721902	ø20				25/450
00721704	-		Qimospike V is a vented spike for vials fitted with Qimo female ♀	Repeated withdrawals from every type of vial	25/300
00721302	-		Qimoharpoon is a closed system transfer device that has an air-intake tube connected to a 50ml balloon filled with sterile air which ensures pressure equalisation	Withdrawal and reconstitution of high-risk medications	10/80
00721301	-				

Syringe


Vygon code	Image	Description	Indication	Quantity box/case
00721091		Qimo male ♂ is a disinfected closed needle-free safety connector. It is a safety captive adaptor for male Luer-lock hubs. Once attached, it cannot be disconnected. It connects only to Qimo female ♀	To secure syringes used for medication preparation	50/900

TRANSFER DEVICES FOR MEDICATIONS/SOLUTIONS



IV fluid bag extension sets

Vygon code	Image	Description	Indication	Quantity box/case
00722101		IV fluid bag extension set with a lateral injection site fitted with Qimo female ♀. Qimo male ♂ secures the distal end of the extension tube	Dilution and administration of all kind of medications	25/150
007221011		IV fluid bag extension set with 0.2µm filter and lateral injection site fitted with Qimo female ♀ . Qimo male ♂ secures the distal end of the extension tube	Dilution and administration of solutions which need to be filtered (e.g. Paclitaxel)	25/150






Female Luer-lock connector

Vygon code	Image	Description	Indication	Quantity box/case
00722001		Female Luer-lock to female Luer-lock connector; with Qimo female ♀ on both ends of the connector.	To transfer medications from syringe to syringe.	50/600


Elastomeric pumps

Vygon code	Image	Description	Indication	Quantity box/case
Option 1 00721004		Set with Qimoprime, Qimo male ♂ and Qimo female ♀ in the same packaging	Can be used to secure the injection site of the elastomeric pump and the distal end of its tubing, but other applications are possible	50/600
Option 2	00721002	Qimo female ♀ is a multipurpose, closed, needleless and disinfectable safety connector; supplied in individual packaging	Qimo female ♀ can be used to connect the tubing of the distal end of the elastomeric pump secured with Qimo male ♂, to access an implantable port or any other patient vascular device	50/1200
	00721001		Set supplied with Qimo male ♂ and Qimoprime in the same packaging	Can be used to secure distal end of elastomeric pump tubing. Need to be connected to a Qimo female ♀


INFUSION DEVICES

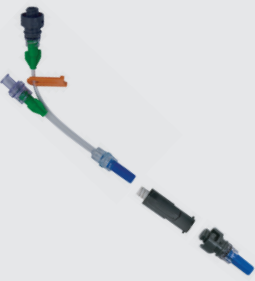Vygon code	Image	Description	Indication	Quantity box/case
007227101		Qimopump is a safe adaptor for pump infusion sets. This device is fitted with 1, 2 or 4 Qimo female ♀ connectors, a spike with an air-venting port (with a cap) and a universal tube to connect to dedicated pump sets	To secure the dedicated pump infusion sets and connect the bag containing medication in a safe manner	25/150
007227201				
007227401				
00722102		IV fluid bag extension tube with a lateral injection site fitted with Qimo female ♀ and an air-venting port (with a cap)	Preparation and administration of non-toxic drugs (premedication such as anti-emetics, corticoids, etc.)	25/150
00085502		Tubing set is an extension tube with vented spike and a male Luer-lock connector with a non-return valve and a priming cap	Extension tube for drug administration	25/150
00721001		Set supplied with Qimo male ♂ and Qimoprime in the same packaging	To secure both the male Luer-lock and the female Luer-lock end of infusion devices	50/600
007231221		Qimooctopus is a safe double-lumen IV extension tube allowing continuous infusion, intermittent injection and blood sampling. Its lateral line is fitted with Qimo female ♀ Qimo male ♂ has to be attached to the distal end of Qimooctopus. Qimoprime allows the priming of the line once Qimo male ♂ has been attached	Safe emergency access. Easy direct injection with a syringe close to the patient	50/300

Vesical instillation

Vygon code	Image	Description	Indication	Quantity box/case
00721006		qimoadapt is a safe adaptor for the connection of a syringe on the funnel of a tube. Qimo female ♀ allows to connect a syringe secured with Qimo male ♂	Safe bladder instillation	50/600

VASCULAR ACCESS FOR PATIENTS

Vygon code	Image	Description	Quantity box/case
007242720		Qimoperfusafe is a safety infusion set with Huber needle fitted with a Qimo female ♀ for the access of implantable ports	10/90
007242725			
007242920			
007242925			
00721002		Qimo female ♀ is a multi-purpose, closed, needle-free and disinfectable safety connector	100/1800

Implement Qimono today

Contact us today to see how Qimono can help your Oncology Team prepare, transfer and infuse chemotherapy drugs safely.

Call **01793 748830** or email marketing@vygon.co.uk

References

1. Simons A, Toland S (2015) Closed systems for drug delivery; a necessity, not an option British Journal of Nursing, (IVTherapy Supplement), Vol 24, No 14
2. Safe Handling of Cytotoxic Drugs in the Workplace, HSE 2015
3. Simons A, Toland S (2017) Perceived adverse effects from handling systemic anti-cancer therapy agents British Journal of Nursing, (IVTherapy Supplement), Vol 26, No 16

Further reading

1. Preventing Occupational Exposures to Antineoplastic and other Hazardous Drugs in Health Care Settings, National Institute for Occupational Safety and Health (NIOSH), September 2004
2. ISOPP Standards of Practice: safe handling of cytotoxics - Section 7 Medical devices, International Society of Oncology Pharmacy Practitioners (ISOPP), J Oncol Pharm Pract, 2007.
3. Infusion Nurses Standards of Practice, Page S32, section 27, Practice Criteria A, The Infusion Nurses Society (INS), 2011.
4. Choosing the best design for Intravenous needle-free connectors to prevent bloodstream infections, William R, Jarvis, MD, Infection Control Today, 28th July 2010.
5. Combattre les TMS par l'ergonomie, Florence Martinache, Industrie Pharma magazine. 1st May 2012.
6. Troubles Musculo-squelettiques : analyse critique des équipements hospitaliers, I. Juret, V. Josselin, S. Moisan, Chinon (CHU Angers, France), Section 3.2, p.45.
7. Efficacy of the Valve Systems of Needle-Free Closed Connectors, report 67-08-A, The Health Protection Agency UK, 21 May 2009.
8. An evaluation of Bionector Microbial Integrity, report 65-07, The Health Protection Agency UK, 28th November 2007.
9. Validation of the Manual Disinfection Procedure after multiple handling, NAMSA, report 178835, 28th February 2014.
10. Chemical resistant connection - WP3 QIMO, Vygon SA, December 2015.
11. Qimo male ♀ Luer compatibility - WP2 QIMO, Vygon SA, December 2015.
12. Qimono closed system - WP4 QIMO, Vygon SA, December 2015.

FOR FURTHER INFORMATION, PLEASE CONTACT: info@vygon.co.uk

The specifications shown in this leaflet are for information only and are not, under any circumstances, of a contractual nature. This brochure has been printed on responsibly sourced and sustainable material. To help us reduce our carbon footprint all of our literature is available electronically either from your Product Specialist, on our website or by emailing info@vygon.co.uk

VYGON (UK) LTD, THE PIERRE SIMONET BUILDING, V PARK,
GATEWAY NORTH, LATHAM ROAD, SWINDON, WILTSHIRE SN25 4DL
RECEPTION: +44 (0)1793 748800 WWW.VYGON.CO.UK

 [@vygonuk](https://www.vygon.co.uk)